



StarDrop™

StarDrop 8

# StarDrop™ Worked Example: Generation of New Compound Ideas



Copyright © 2026 Optibrium, Ltd. All Rights Reserved.

Optibrium™, StarDrop™, Inspyra™, Augmented Chemistry®, Cerella™, Glowing Molecule™, WhichP450™, WhichEnzyme™, Auto-Modeller™, Nova™, MPO Explorer™ and Card View® are trademarks of Optibrium Ltd. Matsy™ is a trademark of NextMove Software Ltd., BIOSTER™ is a trademark of iKem Szolgáltató és Kereskedelmi BT, Derek Nexus™ is a trademark of Lhasa Ltd., Surfex eSim3D™ is a trademark of BioPharmics LLC, and SeeSAR™ is a trademark of BioSolveIT GmbH. US Patent Numbers 9,224,098 and 9,367,812

## Contents

Introduction	3
Loading data and model files	3
Generating new compound ideas with Nova	4
Interpreting the results of compound generation	9
Conclusion	13

# Introduction

This example considers a challenge typically faced when a series of promising compounds has been found, but the chemistry needs to be broadened to further improve one or more properties. Focussing on a single property for optimisation often increases the risk of introducing liabilities for other properties, and so in this exercise, we're going to see how StarDrop™ can help to explore new chemistry ideas while minimising this risk.

During this example, we will consider three compounds from a lead series that we would like to try to evolve into a candidate. The compound has a good profile of ADME properties but insufficient inhibition of the target, the Serotonin transporter. In this example, we will use StarDrop's Nova™ module to generate new ideas for compounds to improve the potency while maintaining the balance of other properties.


We will use a model for Serotonin transporter inhibition built with public domain data from the ChEMBL database using StarDrop's Auto-Modeller™ to monitor potency during the exploration.

Nova can generate data sets of many hundreds of thousands of compounds if left to run for multiple generations, and so in this example, we will also explore how this process can be managed to focus on those compounds most likely to achieve our objectives.

## Loading data and model files

- In StarDrop, select **Open** from the **File** menu, find the StarDrop file **Nova lead compounds.sdproj** and select it to open it.


This simple project file contains a data set with three lead compounds for which a range of properties have been predicted. The compounds have been scored using StarDrop's approach to multi-parameter optimisation, called Probabilistic Scoring, against a profile of properties suitable for an orally dosed compound intended for a CNS target.

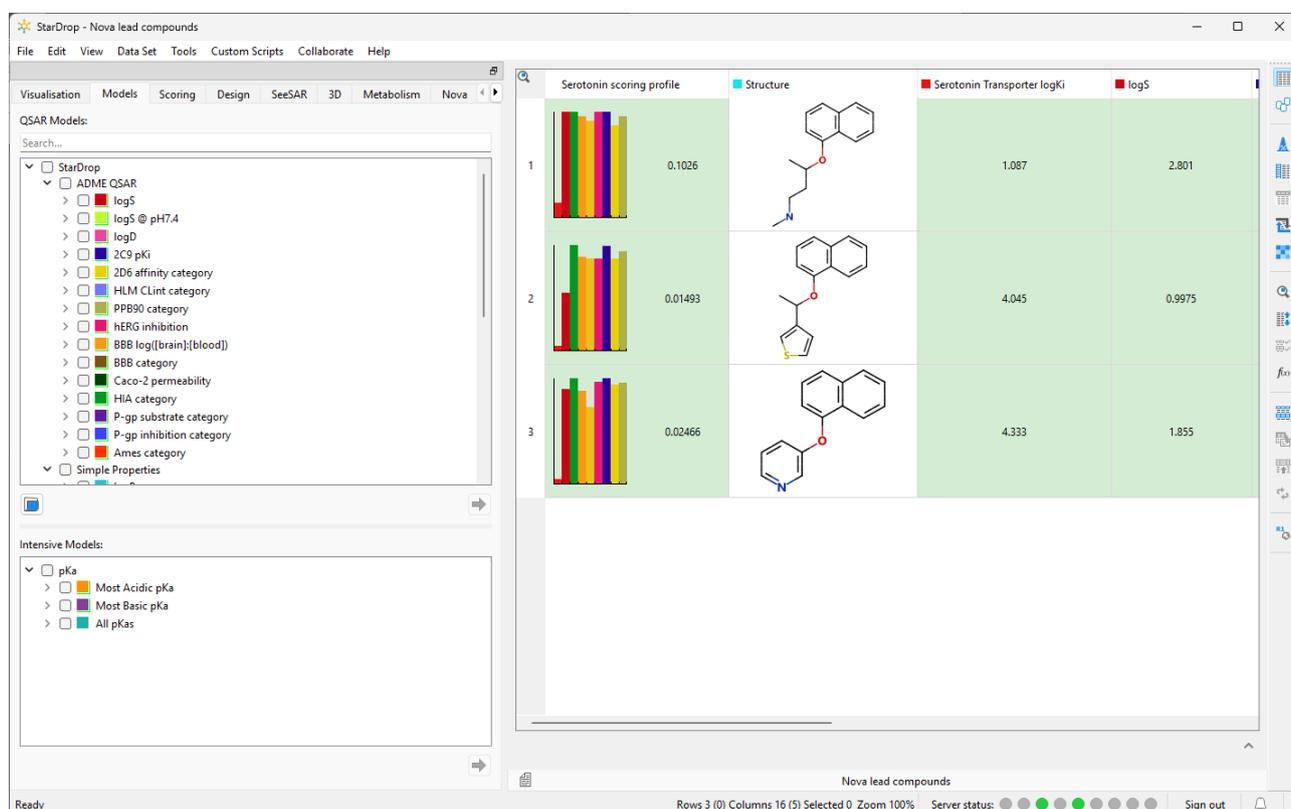

The third column in the data set is called Serotonin Transporter logKi (uM), and the values in this column have been generated by a custom model built using ChEMBL data with the StarDrop Auto-Modeller. These values are used when calculating the overall scores; therefore, to calculate this property and scores for any new compounds, we need to load this model into StarDrop.

- Click the **Open** button  at the bottom of the **QSAR Models** area and open the model file **Serotonin\_Transporter\_logKi.aim**.

You will now see that a new model **Serotonin Transporter logKi** is displayed in the **Custom** branch in the **QSAR Models** area. StarDrop will now be able to use this model in the prioritisation of compounds during the Nova idea generation process.

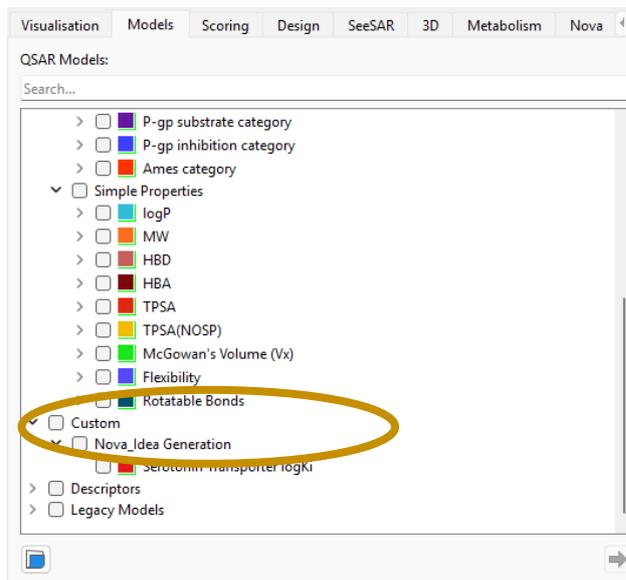

## Generating new compound ideas with Nova


- Click on the **Nova** tab and select all three rows in the data set by clicking in the top left corner of the data set.
- Click the  button to start the Nova wizard.

Row	Serotonin scoring profile	Structure	Serotonin Transporter logKi	logS
1	0.1026	<chem>CC(C)OC1=CC=C2C=CC=C12</chem>	1.087	2.801
2	0.01493	<chem>CC(C)OC1=CC=C2C=CC=C12</chem>	4.045	0.9975
3	0.02466	<chem>C1=CC=C2C=CC=C12</chem>	4.333	1.855

The Nova dialog will be displayed. On the first page, you can choose whether to apply Chemistry Transformations, enumerate a library or generate Matched Series Suggestions.

- For this example, choose **Chemistry Transformations** and click the **Next** button.



**Note:** Additional worked examples are available for the other approaches.

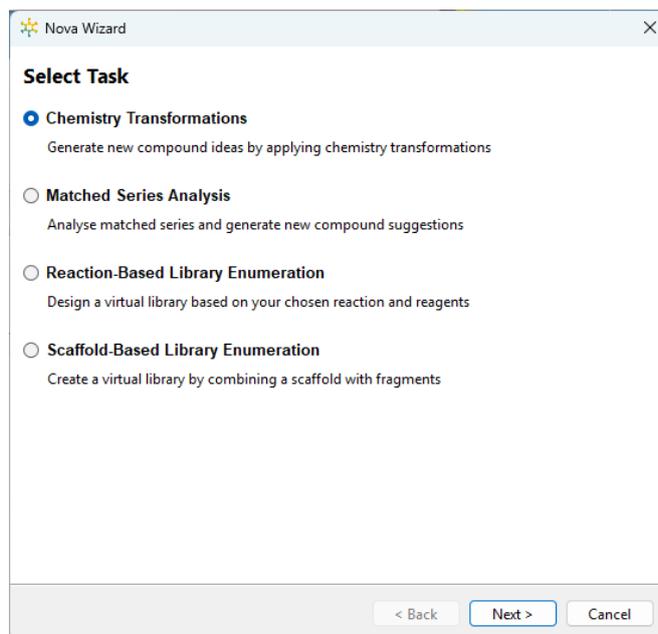

The **Specify Input Structure** page of the dialog displays the structures we selected from the data set.

If you wish to indicate regions of the molecules that must remain unchanged, then you can draw around them with the mouse to create a mask.

Selecting the individual molecules on the left enables you to define individual masks for each molecule. A new compound will be generated and displayed.

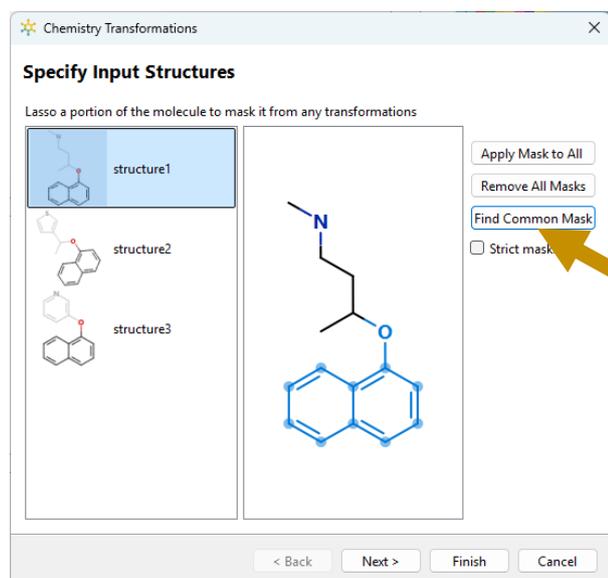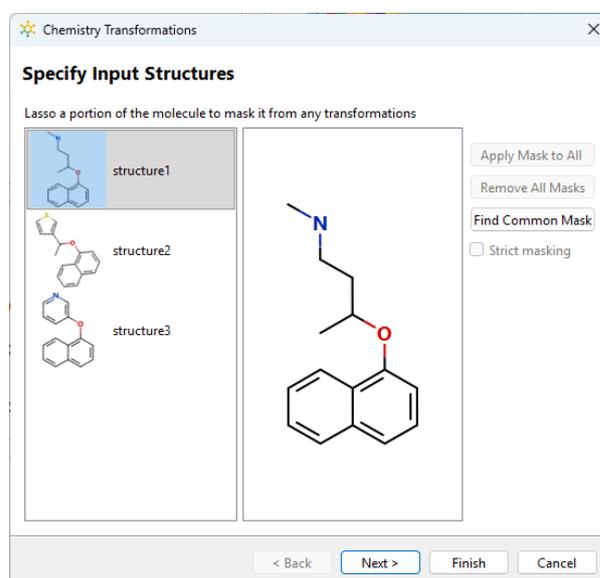

In this example, all the molecules contain a common sub-structure.

- Rather than drawing this for each compound, click the **Find Common Mask** button.




The **naphthol** group will now be masked on each compound.

- Click the **Next** button to go to the **Select Transformations** page.

On this page, we can select which transformations we want to apply. They are grouped to make it easy to pick certain types of transformation. The **Manage** button enables you to reorganise the transformations or import your own.

**Note:** Nova is showing the entire library of transformations, including those that do not apply to our specific structures. If you click the **Show applicable** option, Nova will search the database of available transformations to find those that apply to the input compounds. The **Local** transformations are those provided with the Nova module. If you have a license for the **BIOSTER™** module, you will also see a database of >29,000 precedented transformations manually curated from the literature.

You can open the groups of transformations by clicking on the arrow next to the group to browse the individual transformations, as shown to the right.

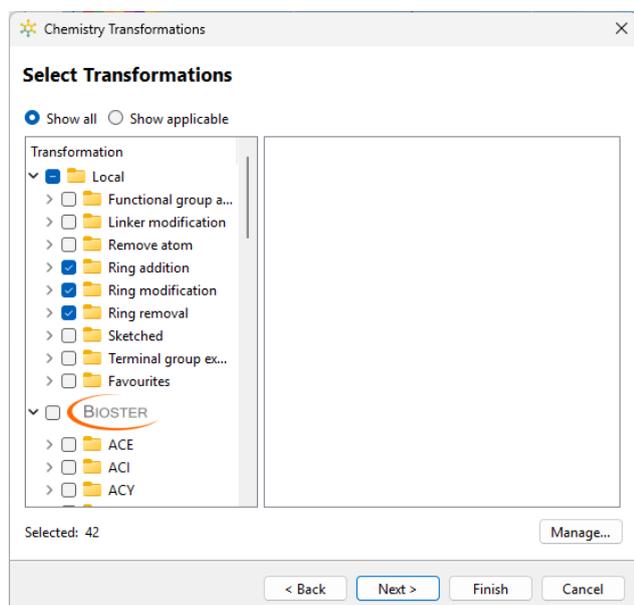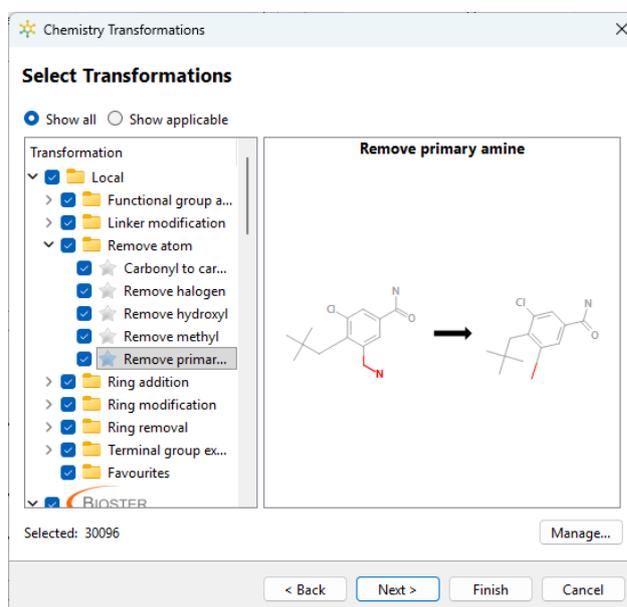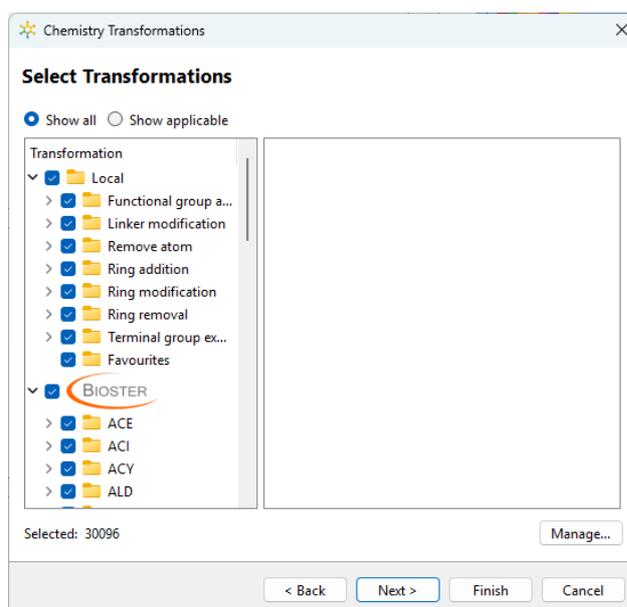

- In this quick example, we're only going to select the three transformation groups relating to rings: **Ring addition**, **Ring modification** and **Ring removal**. Untick the boxes next to the other groups and the **BIOSTER** library, as shown to the left.

- With **only** these three items selected, click the **Next** button to go to the **Control Output** page.


Here, we can choose how many generations of transformations to apply.

- Change the number of **Generations** to **3**.

If we apply multiple generations, then all of the molecules created in one generation will be used as starting points for the next generation. This can result in a very large number of compounds being generated.


- To restrict the number of compounds used as the basis for subsequent generations, tick the **Select subset each generation** checkbox.

First, we choose a **Method**. We can select compounds based on a property, a score, chemical diversity or a combination of these. If a property or score is involved, then we can indicate whether it is desirable for this to be high or low.

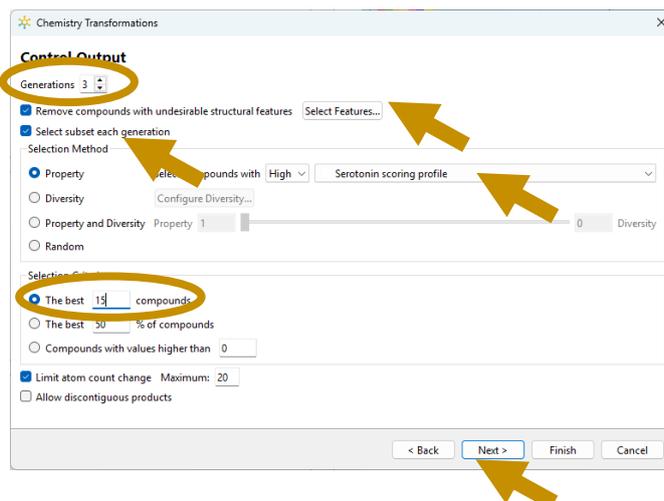

- In this case, we are going to make a selection that is biased entirely towards having a **High** score against the **Serotonin scoring profile**, so select that.
- We can then define the **Selection Criteria**. In this case, choose **The best 15** compounds.
- Lastly, click the **Select Features** button.

On this page, we can indicate structural features we would like to ensure are not present in any new compounds. Compounds that contain any of these will be removed at the end of the process. A number of published sets of filters are available, and by clicking the **Import** button, you can add your own, but in this example, we will just use the **Undesirable fragments** that are already checked.

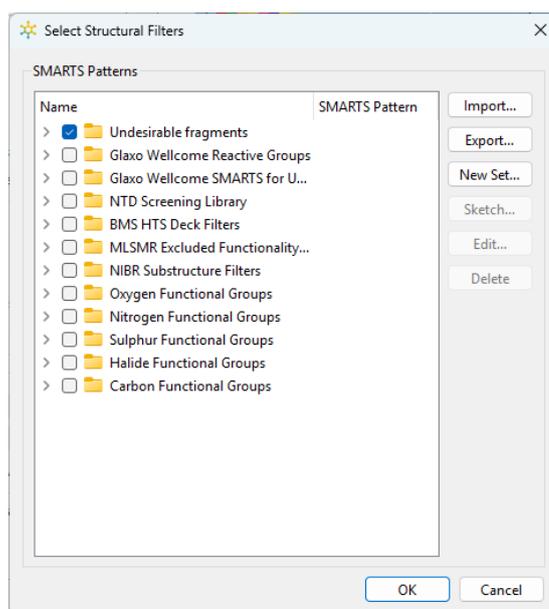


- Click the **Next** button to go to the **Create New Data Set** page.
- On the final page, we can give the Nova session a name, in this case, **Serotonin Example**, and add any other relevant notes. The session name will be used for the new data set that is created.
- Tick the **Show results in Card View** checkbox at the bottom of the page.
- Start the idea generation process by clicking **Finish**.

The Nova process will take a couple of minutes to run. While the process is running, an indicator will show which generation is currently underway, along with the best property or score value achieved so far.


Generation	Structure	Serotonin Transporter logKi	logS
1	<chem>CC(C)OC1=CC=C2C=CC=CC12</chem>	1.087	2.801
2	<chem>CC1=CC=C(C=C1)OC2=CC=CC=C2</chem>	4.045	0.9975
3	<chem>CC1=CC=C(C=C1)OC2=CC=CC=C2</chem>	4.333	1.855

# Interpreting the results of compound generation


Once complete, a new data set called **Serotonin Example** will be displayed in Card View.

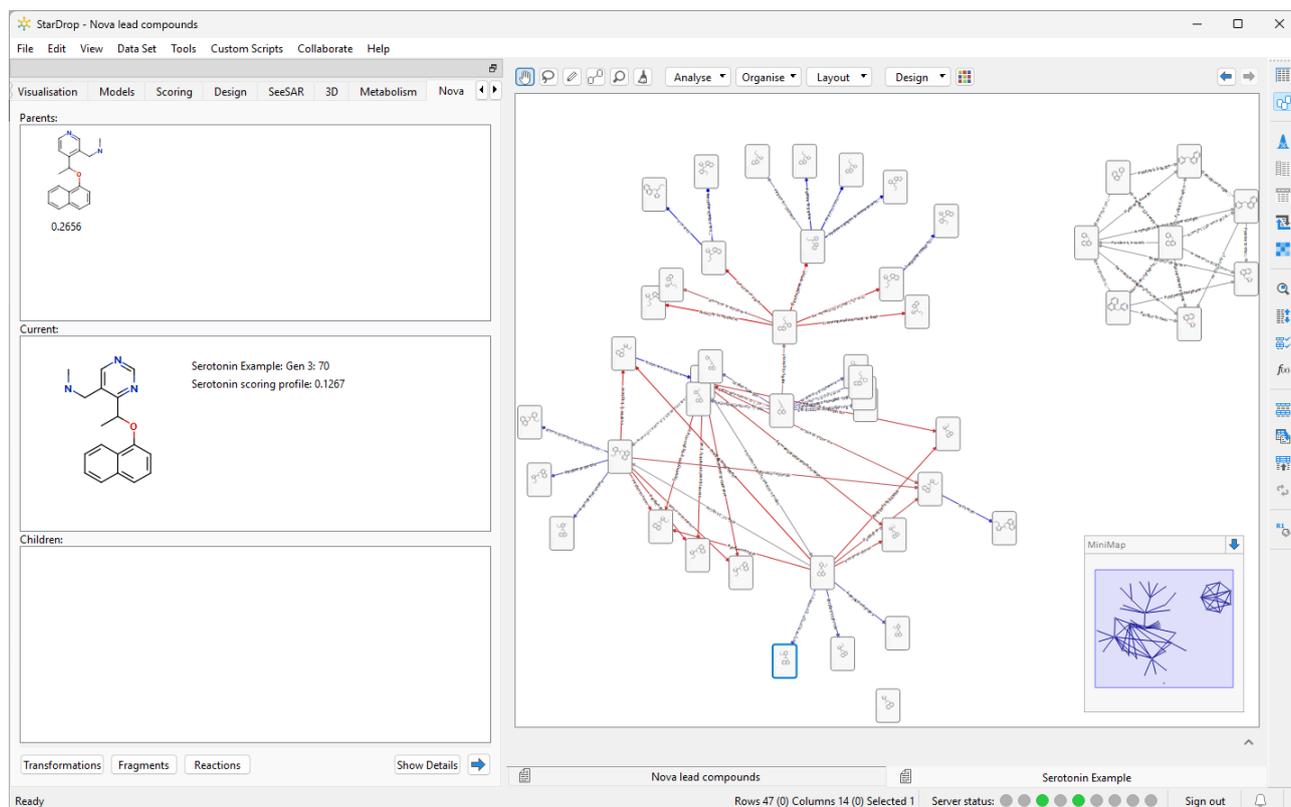


- We can colour the cards to show the different generations of compounds by clicking the **Format** button  at the top of Card View.
- In the **Colour By** dialogue, choose the property **Generation**.

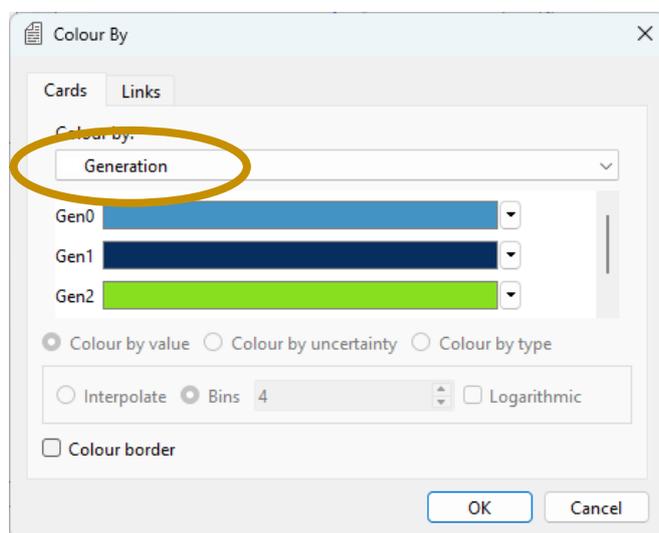


From the resulting display, shown below, we can see that, of the three initial compounds, one (top-right) resulted in the selection of compounds only from the first generation, and none of the subsequent generations had a sufficiently high score to be selected. Another compound (bottom) did not result in any 'children' that were selected, while all the second and third-generation compounds were generated from only one of the initial compounds (centre).

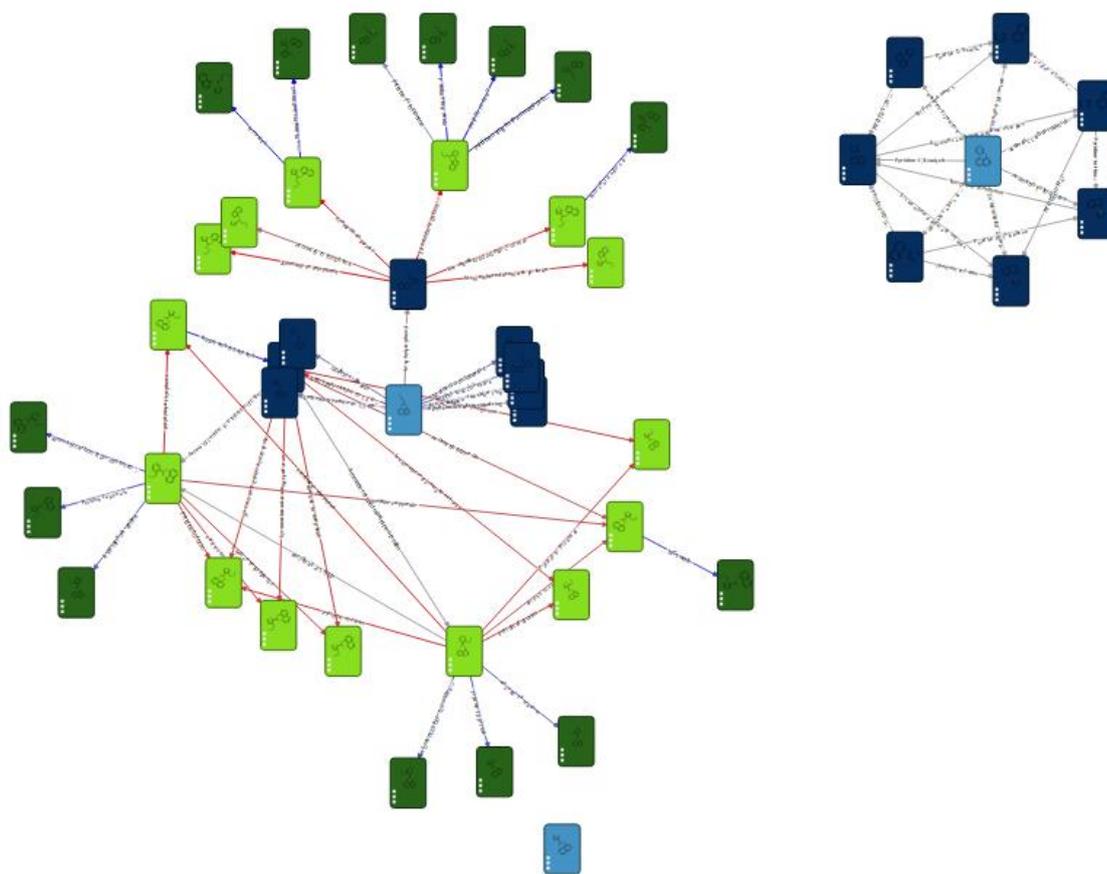

- We can now choose to colour the cards by the overall score by clicking again on the **Format** button  at the top of Card View.
- This time, in the **Colour By** dialogue, choose the property **Serotonin scoring profile**.

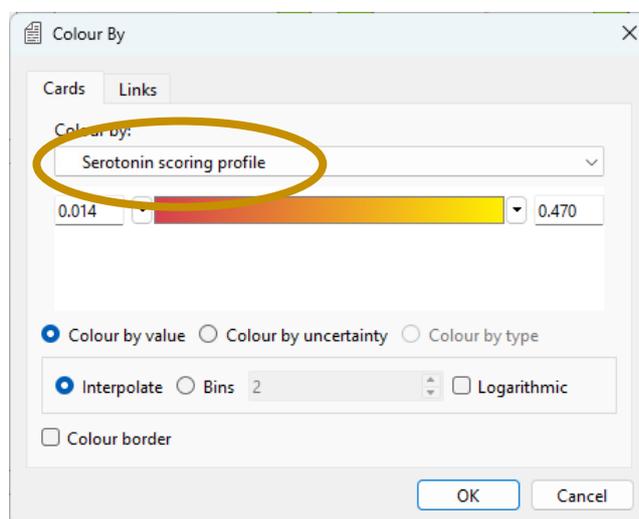


From this display (below), we can see that most of the high scoring compounds are in the second generation (the yellowest cards).

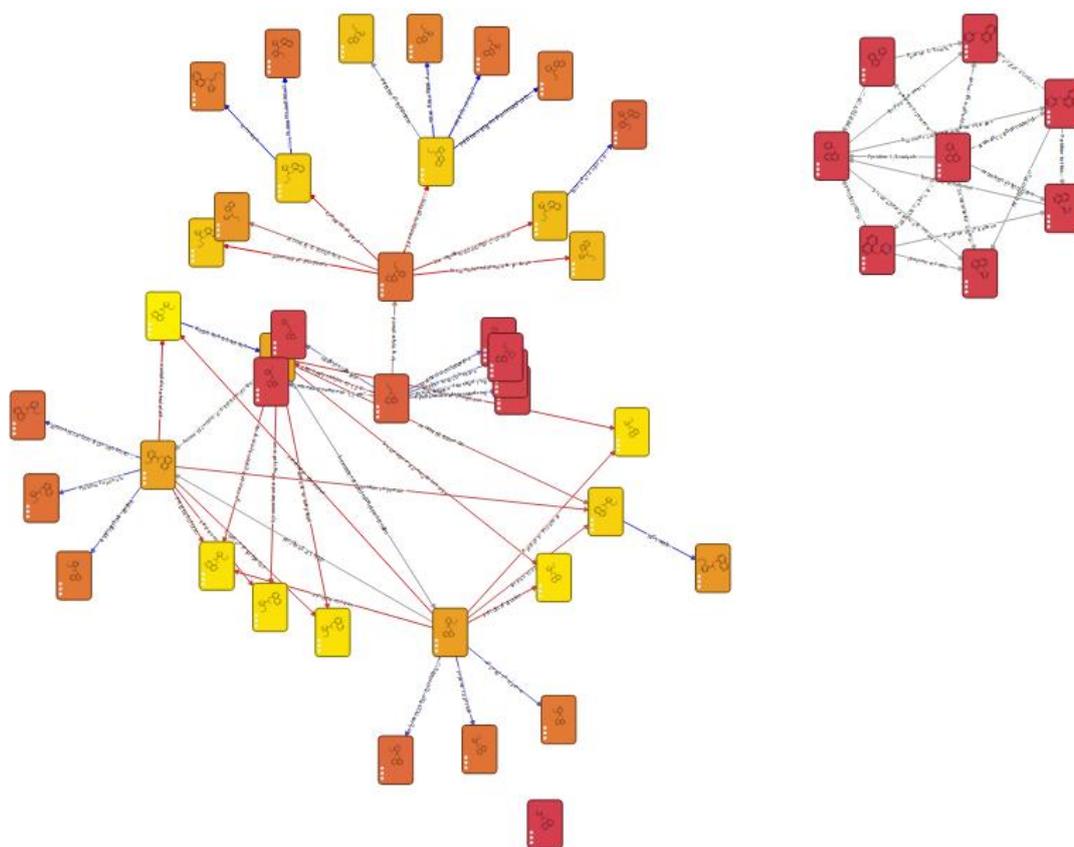

You will also note that the links between cards are coloured by the change in score between the corresponding compounds; a large increase is shown by a red link, a large decrease by a blue link and a small change in grey.

- Zoom in to the compound circled in green above using the mouse wheel (or by simultaneously pressing the **Ctrl** and = keys)

From this, you will note that the biggest improvements in score from this compound (the red links) result from the “Aromatic to thiophene” and “Aromatic to pyrrole” transformations.

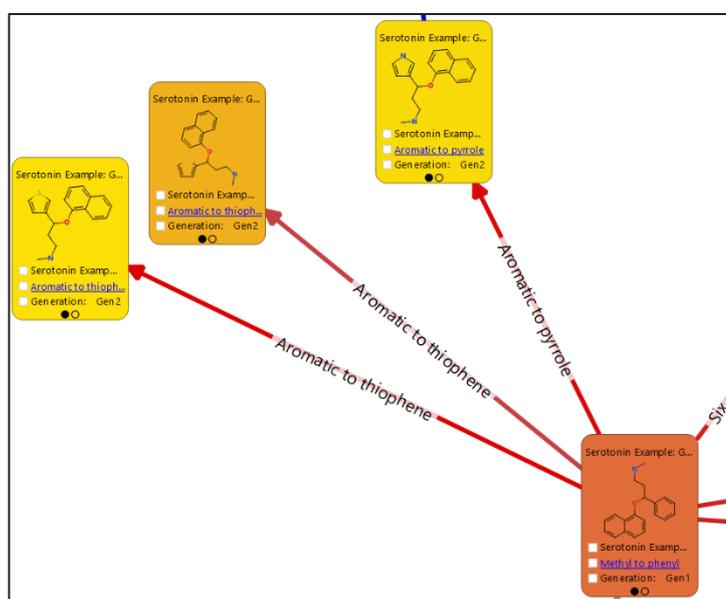

You can see a transformation in more detail by clicking on the hyperlink on the cards. You can also choose to see more or different properties on the cards by clicking on the **Design** menu at the top of Card View.

- Switch to Table View by clicking on the **Table View** button .

The compounds in the data set are sorted to show the compounds with the best property or score at the top. The data set contains a column showing which transformation was applied to generate each compound from its immediate ‘parent’ and a column indicating in which generation this occurred.

- Select the first row.

The screenshot shows the StarDrop interface with a table of generated compounds. The table has the following columns: Structure, ID, Transform, Generation, and Serotonin Transporter Inhibition (gKi). The first row is highlighted in blue, indicating it is selected.

Structure	ID	Transform	Generation	Serotonin Transporter Inhibition (gKi)
	Serotonin Example: Gen 2: 89	Aromatic to thiophene	Gen2	0.2057
	Serotonin Example: Gen 2: 88	Aromatic to thiophene	Gen2	0.357
	Serotonin Example: Gen 2: 90	Aromatic to thiophene	Gen2	0.357
	Serotonin Example: Gen 2: 81	Six membered aromatic to ...	Gen2	0.4585
	Serotonin Example: Gen 2: 80	Six membered aromatic to ...	Gen2	0.4585
	Serotonin Example: Gen 2: 82	Six membered aromatic to ...	Gen2	0.4585
	Serotonin Example: Gen 2: 86	Aromatic to pyrrole	Gen2	0.3068
	Serotonin Example: Gen 2: 33	Aromatic to pyrrole	Gen2	0.3811
	Serotonin Example: Gen 2: 30	Six membered aromatic to ...	Gen2	0.4833
	Serotonin Example: Gen 2: 35	Aromatic to thiophene	Gen2	0.3923

The left panel shows the 'Parents' section with three chemical structures and their scores: 0.281, 0.2715, and 0.2656. The 'Current' section shows the selected compound (Serotonin Example: Gen 2: 89) with a Serotonin scoring profile of 0.4694. The 'Children' section is currently empty.

The first row contains the top-scoring compound, and we can see this compound's parent molecules and any subsequent children in the viewer on the left. This compound was created in generation 2.

Sometimes, the same molecule can be generated from different parent molecules. In this case, all three parent compounds could be transformed with the "Aromatic to thiophene" transformation to produce this molecule. Notably, all three parent compounds have scores of around 0.27 or higher, but the introduction of the thiophene ring significantly improves this.

**Note:** Structures that are greyed out have been removed from the dataset either because they are duplicates of a molecule previously invented or because they contained one of the structural filters in the final filtering step.

Not only have we identified new compound ideas with a higher likelihood of activity at the target receptor, but we have also identified compounds that are medically relevant. Notice that row 17 is actually the drug Duloxetine, a serotonin-norepinephrine dual reuptake inhibitor. Its isomer at the 3-position of the thiophene ring, CHEMBL142028, which has also been shown to have similar serotonin reuptake activity, can be found in row 10.

StarDrop - Nova lead compounds

File Edit View Data Set Tools Custom Scripts Collaborate Help

Visualisation Models Scoring Design SeeSAR 3D Metabolism Nova

Parents: 0.1436

Current: Serotonin Example: Gen 2: 34  
Serotonin scoring profile: 0.2443

Children:

Transformations Fragments Reactions Show Details

Structure	ID	Transform	Generation	Serotonin Transporter I <sub>0</sub> gK <sub>i</sub>
	Serotonin Example: Gen 2: 35	Aromatic to thiophene	Gen2	0.3923
	Serotonin Example: Gen 3: 22	Pyridine to thiophene	Gen3	0.5399
	Serotonin Example: Gen 2: 28	Six membered aromatic to ...	Gen2	0.5225
	Serotonin Example: Gen 2: 27	Six membered aromatic to ...	Gen2	0.5225
	Serotonin Example: Gen 1: 7	Secondary carbon to 1,2-p...	Gen1	0.5607
	Serotonin Example: Gen 2: 77	Aromatic carbon to pyridin...	Gen2	0.6724
	Serotonin Example: Gen 2: 78	Aromatic carbon to pyridin...	Gen2	0.6726
	Serotonin Example: Gen 2: 34	Aromatic to thiophene	Gen2	0.5463
	Serotonin Example: Gen 3: 98	NC-switch	Gen3	0.7397

Ready Nova lead compounds Serotonin Example Rows 47 (0) Columns 25 (11) Selected 1 Zoom 100% Server status: Sign out

## Conclusion

Further details about using Nova in conjunction with StarDrop's Probabilistic Scoring can be found in the following publication "Applying Medicinal Chemistry Transformations and Multi-parameter Optimization to Guide the Search for High-Quality Leads and Candidates", Segall et al., J. Chem. Inf. Model., 2011, 51 (11), pp 2967–2976, also available as a preprint from:

<https://optibrium.com/knowledge-base/applying-med-chem-transformations-and-multi-parameter-optimisation/>

This example does not explore the potential of applying over 29,000 precedented transformations using the BIOSTER module. A short introductory video to the BIOSTER module can be found here:

<https://optibrium.com/products/stardrop/modules/bioster/>

If you have any questions, please feel free to contact [support@optibrium.com](mailto:support@optibrium.com).